

MODEL-med
international pty. ltd.


PPH KIT

ABN 23 492 360 470

1296 Glenhuntly Rd, Carnegie 3163
Melbourne Australia

Phone +61 (0)418 506 308 www.modelmed.com.au

Sometimes the rubber in Model-med models will be dented or flattened if left for a long time in one place.

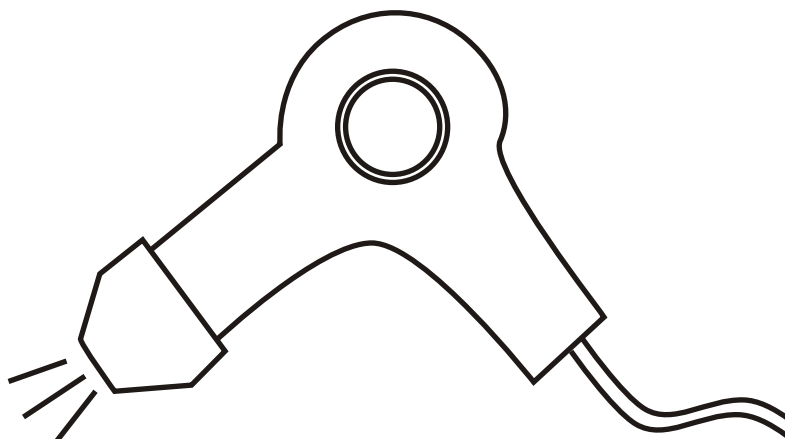

If you find that the rubber is deformed, gently warm the rubber with a hairdrier. But be careful not to overheat the rubber. A gentle heat is all that is needed.



This will bring the rubber back to its original shape.

MODEL-med
international pty. ltd.

ABN 23 492 360 470

1296 Glenhuntly Rd, Carnegie 3163
Melbourne Australia

Phone +61 (0)418 506 308

www.modelmed.com.au

Post-Partum Haemorrhage Kit Care and Handling

The Model-med PPH Kit is simple to use.

There are two tubes coming out of the uterus. One is for a low bleed and the other is for a higher bleed in the uterus. Choose which tube you wish to use.

Now connect the tube to the pump as follows:-

Place the hose clamp over the tube coming from the nozzle of the pump.



Now insert the tube from the uterus all the way into the tube from the pump and gently tighten the hose clamp to prevent artificial blood from leaking out.



Make sure the birth canal of Sophie's Mum is clean and dry. Place the uterus inside the birth canal until the opening of the uterus is at the opening of the vagina. Rather than clamping the Sophie's Mum base to the table, clamp the clamps to the table and butt the base up against them. Then raise the back of the model by approximately 150mm and rest it on a support such as a folded towel. This will assist the blood to flow out of the vagina rather than pooling inside the uterus.



Unscrew the top of the pump and fill the bottle with Model-med Artificial Blood to a maximum of 2 litres. Tighten the top back onto the bottle and pump approximately 80 times. This pressure is required because the blood is very viscous. Depress the green trigger to make the blood flow. The small grey tab above the trigger is a lock for the trigger. Push it forward to lock the trigger open and pull it back towards you to release it.

CLEAN UP

Model-med Artificial Blood is non-toxic and will wash out. You may find it will stain your hands red but this disappears quite rapidly. When cleaning the Sophie's Mum you will find that the red stain requires more cleaning than usual. You can use Kerosene or a mild citrus-based cleaner. Citrus-based cleaners come in different strengths, so test it first on a small patch of skin. In Australia we recommend Citro Clean from Chem Pack if you can get it. Wash out the pump and the uterus with warm water.

As usual you should clean the uterus with detergent and warm water, pat it dry, and then powder it with talcum powder to help preserve the rubber.

MODEL-med
international pty. ltd.

ABN 23 492 360 470

Melbourne Australia

Phone +61 (0)418 506 308

Model-med International Pty Ltd Warranty

Our products are made from a soft, yet tough, material. They will last if they are treated like a person. That is, gently.

Model-med International Pty Ltd offers a one year warranty to the original purchaser from the date of purchase on our models when they are handled, operated and maintained in accordance with the printed instructions accompanying them. This Limited Warranty shall be limited to the repair or replacement of parts that prove defective under normal use and service and which on examination shall indicate, to Model-med's satisfaction, they are defective. If Model-med confirms the defect and approves the claim, Model-med will elect to repair or replace, or have parts repaired or replaced without charge by a nominated representative of Model-med 's choosing. We may ask for photographic evidence of the damage to the model. The Warranty does not cover any costs of shipping. Shipping costs will be paid by the Customer.

This Limited Warranty does not cover any failures or operating difficulties due to accident, abuse, misuse, alteration, misapplication, vandalism, improper maintenance, or failure to perform normal and routine maintenance.

We will repair a model if it fails within one year of purchase providing it has not been mistreated or abused. We define mistreated or abused in, but not limited to, the following ways:-

If the model is cut, punctured or sliced by any object, including but not limited to forceps, knives, or any other implements.

Damage caused by fingernails is not covered by this warranty.

If any lubricant other than a water-based lubricant is used on the model. Only recognised water-based surgical lubricants should be used.

If any chemicals other than those described in the care and handling documents are used on the skin, or if the chemicals described in the care and handling documents are used in quantities greater than those recommended.

Hot objects should not be placed on or around the model. We define hot as being above 50 degrees celcius.

If you believe you have a claim under warranty you must first contact your local distributor for instructions on how to proceed with the claim.

Model-med International Pty Ltd shall not be liable under this or any implied warranty for incidental or consequential damages.